


Alcovisor Jupiter

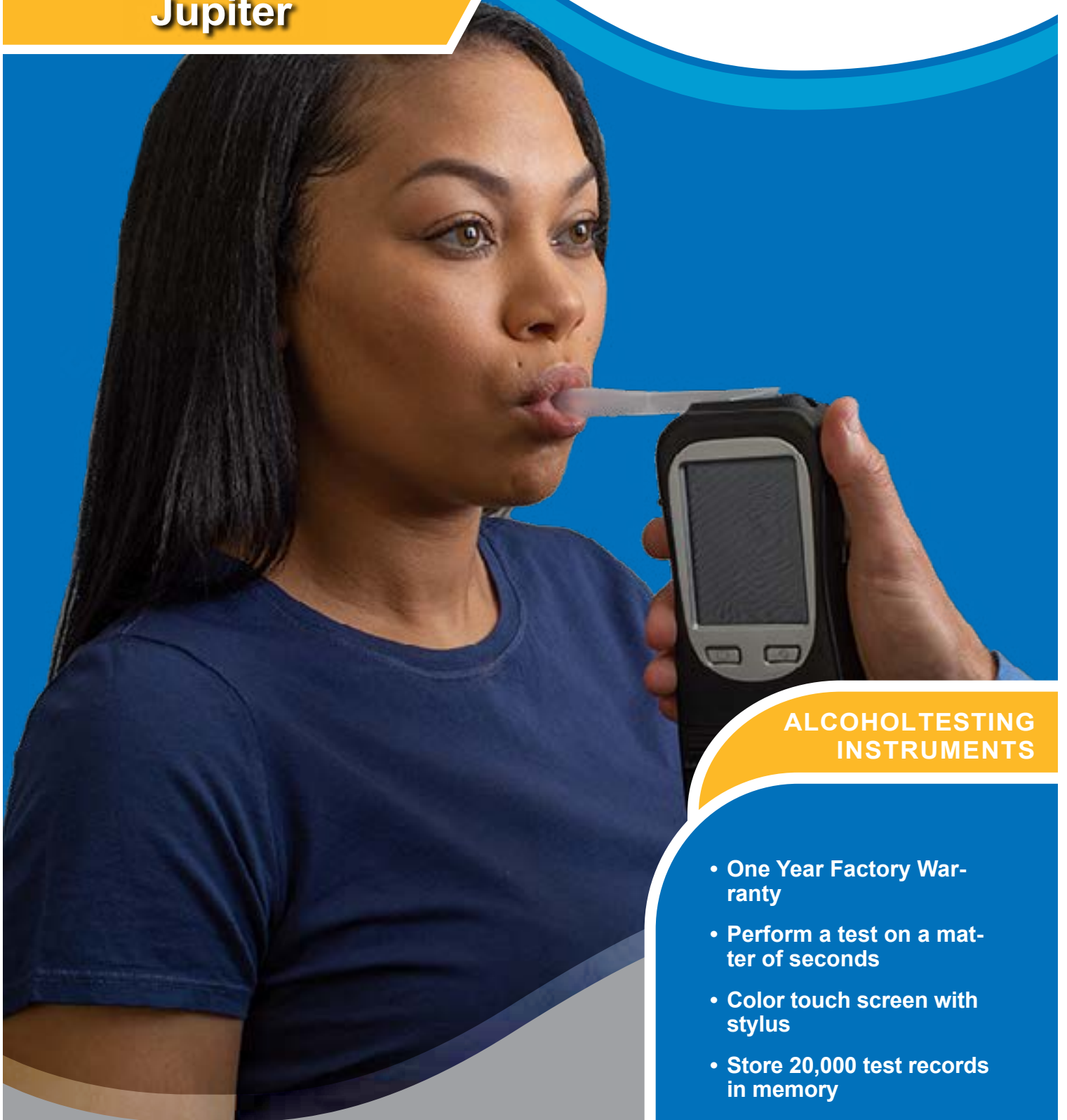

ALCOHOL TESTING INSTRUMENTS

- One Year Factory Warranty
- Perform a test on a matter of seconds
- Color touch screen with stylus
- Store 20,000 test records in memory
- Built-in thermal printer
- NHTSA approved Evidential Breath Testing Device

800.742.4909
kahntactmedical.com

Alcovisor Jupiter Breathalyzer



Quick DOT Testing

The Jupiter system is approved by NHTSA as an Evidential Breath Tester and approved for DOT workplace testing. Fast and easy to use, the Jupiter conducts a breath alcohol test in less than 60 seconds.

Jupiter Product Benefits & Features

- Able to store 20,000 test records in its memory
- Fast testing in a matter of seconds using the Quick Test mode
- :15 countdown timer can be enabled to count down the :15 wait period
- Test with data input mode allows you to define data fields you wish to capture and print
- Automatic and manual breath sampling enabled
- Color LCD touch screen with stylus
- Evidential accuracy
- Built-in thermal printer
- Audible beeping as subject provides a breath sample let's you know air is flowing through the device
- Visual bar on the LCD screen shows you how much breath the subject is providing
- Automatic Air Blank taken at the start of each test
- USB PC Connection available with optional upload software allows you to upload test records to a computer

Stores Thousands of Test Results

The robust memory of the Alcovisor Jupiter allows the instrument to store up to 20,000 test results. The operator can browse the test results stored in the instrument and re-print selected tests, or upload results to a computer with the optional upload software.

RELATED ITEMS AND SUPPLIES

Item	Description	Price
388-JUP	Mouthpiece for the Jupiter and Mercury system	\$.41 @ 1-499 \$.39 @ 500-999 \$.36 @ 1,000+
388-THERM	Alcovisor Printer Paper (4 rolls)	\$13.50
324-KT	Alcohol Gas Tank, 105 liter	\$170.00
382-04KT	Alcohol Gas Tank, 34 liters	\$123.00
328	Regulator	\$112.00
371	DOT Alcohol Test Forms	\$0.25 ea.
372	Non-DOT Alcohol Test Forms	\$0.25 ea.
342-KT	Yellow Tamper Tape (1 box)	\$25.00

STANDARD PRODUCT

Item	Description	Price
201	Jupiter System Includes: Device with printer, standard carrying case, 105 mouthpieces, batteries, printer paper, passive sampling cup, manuals, 1-yr warranty	\$1,515.00
201-KIT	Jupiter System (DOT Package) Includes: Device with built-in printer, large carrying case, 105 mouthpieces, batteries, printer paper, passive sampling cups, manuals, 34-liter gas tank, 100 DOT and non-DOT alcohol testing forms, box of tamper tape, 1-yr warranty	\$1,810.00

SPECIFICATIONS

Sensor type: Electrochemical Oxidation (fuel cell)

Range of measurement: .000 - .400 BrAC

Accuracy and precision: U.S. DOT approved for evidential use

Breath sampling: Automatic with manual override

Response time: 3-5 seconds negative results; 5-10 seconds positive results

Power requirement: 7.4V rechargeable Lithium battery (included)

Computer capability: Optional upload software

Calibration: Microprocessor controlled, using gas or wet bath simulator

Memory: 20,000 tests

Print capability: With built-in thermal printer

Working temperature range: -10 to 50°C

Warranty: One year parts and labor

For Information & Orders Call: **800.742.4909**
Order 24 Hours A Day At: www.kahntactmedical.com
Fax Orders To: **708.449.7832**

Quality Medical Devices Since 1987
Knoxville, TN 37939